



# Freddo

SKIN COOLING DEVICE

A PATIENT WITH NO PAIN  
IS A SATISFIED PATIENT



# WELCOME TO FABINJECT LLC

We are a U.S. based company focused on serving the medical industry. We create innovative products with advanced technologies, committed to maximum excellence, aiming to exceed the needs of patients and healthcare professionals. **QUALITY, ACCESSIBILITY AND RELIABILITY** are the pillars of our company since its foundation in April 2023, our main focus is to facilitate and streamline the work of professionals in the areas of medicine, health, aesthetics and beauty, thus improving people's lives.



## ▶ INNOVATIVE DNA:

Fabinject LLC brings into the North American market the DNA of Fabinject Technology, a company with more than 30 years of experience and innovation in the Brazilian and South American markets, with award-winning and acclaimed products by the sector. Always attentive to market demands and focused on constant research, we are at the forefront of the medical device industry, always innovating through state-of-the-art technologies and clinically proven methods.

Our products are recognized with superior performance, safety and efficacy.

We are proud to contribute to the well-being of our patients through our ability to develop innovative and effective solutions.

## ▶ OUR COMMITMENT:

We value trustful relationships.

We are committed to building a successful, solid and lasting partnership with you.

Our goal is to support you and meet your needs and those of your patients, always providing the best experience for the successful achievement of your results.



## ▶ PROVEN EXCELLENCE:

With an office and manufacturing facility located in Winter Park, Florida, FABINJECT LLC stands out for producing high-performance medical devices that offer maximum **QUALITY, RELIABILITY, AND DURABILITY**.

Our passion and responsibility is for the success of our customers and the care they provide to patients through our products. We always seek to improve the work of health professionals and the quality of life of the people who use our devices.

## ▶ FOCUS ON RESULTS:

Fabinject LLC is committed to helping you improve your bottom line by providing healthcare technology products that are affordable, high-quality, and reliable.

For more information about our products and services, please contact us or schedule a visit at our factory.



### THE ULTIMATE PLUG & PLAY

Patented automatic drying system!

### 0,56 kW

Max. power consumption during therapy!

### READY TO USE

In 5 minutes!

### NO CONSUMABLES

Less costs!

### H 28" / W 19" / L 20"

Housing dimensions.

### 120 lbs (55 kg)

Weight.

### PROMPT DELIVERY

We work with reliable carriers!



\*Optional articulating arm.

# THE ULTIMATE PLUG & PLAY !



## UNCOMPARABLE EFFICIENCY WITH LOWER ENERGY CONSUMPTION



Freddo offers the same cooling power with less energy consumption. Designed for maximum efficiency, the equipment produces cold air up to -35°C through the hose, ensuring maximum productivity while minimizing energy consumption, effectively doing more with less energy used.



## PATENTED AUTOMATIC DRYING SYSTEM



Unlike other skin cooling machines, Freddo does not require to daily empty the condensed water, nor does it have any leaks or liquid build-up. This unique patented system offers exclusive convenience and efficiency.



## NO CONSUMABLES - LESS COSTS AND MORE PRACTICALITY



Freddo is a cryotherapy system without consumables such as condensation collectors or filters. Designed as a plug & play solution, it is easy to operate and repair, eliminating recurring costs and bringing greater savings over time.



## PROMPT DELIVERY - THE AGILITY YOU NEED



We know that when purchasing equipment, professionals want to start using it without delay. Our equipments and spare parts are available for immediate delivery, whenever necessary.





## USER FRIENDLY CONTROL PANEL

We developed an intuitive panel with only four buttons: On/Off, Play/Pause, Timer and Airflow.

This means fewer distractions and more time dedicated to patient treatment, whether adjusting the parameters of a laser device or performing procedures such as botulinum toxin application and fillers.

### Orange flat tip:

Especially useful for sclerotherapy.



### White fine tip:

Useful for working with laser or IPL (does not absorb light energy).




## READINESS AND EFFICIENCY

With a simplified panel and a few keys, Freddo is ready to use in just 5 minutes after turning it on.

This agility allows professionals to focus on their patients, not on the equipment.



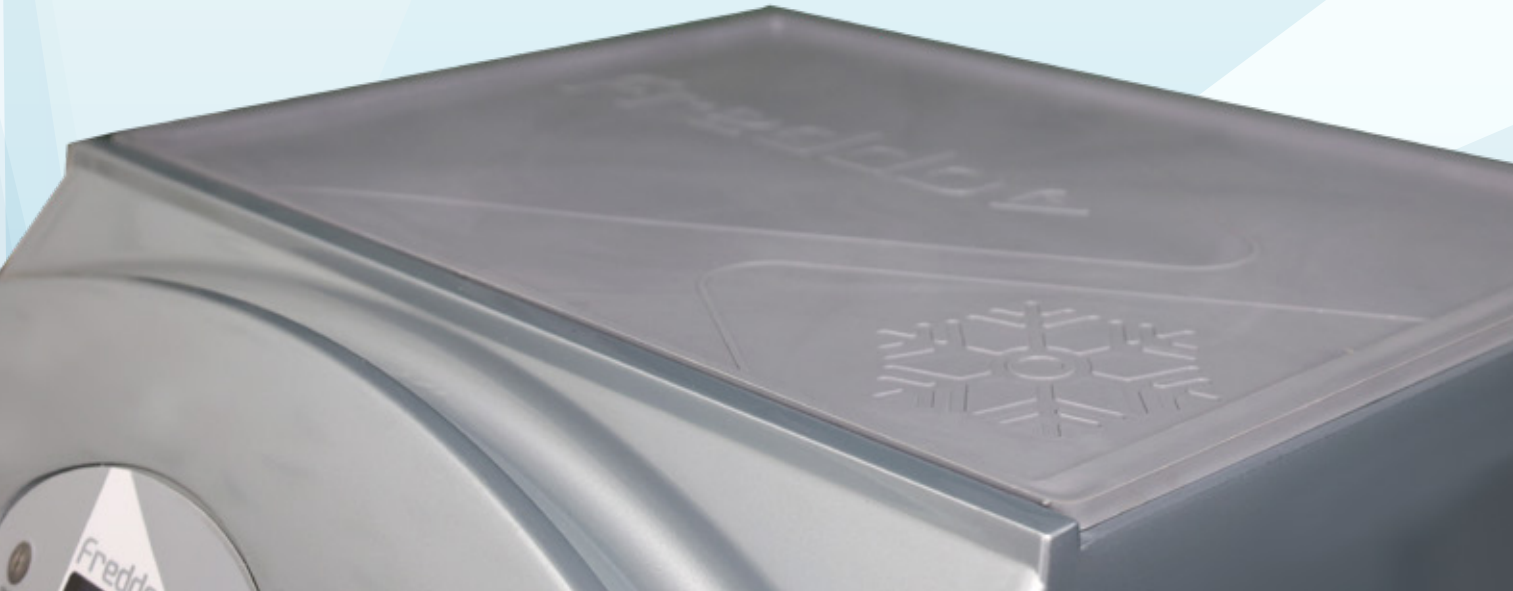




### NO GLASS TOP: MORE SAFETY AND PRACTICALITY

Glass breaks! Our platform is robust and allows us to accommodate other therapy devices and small accessories, without unnecessary risks.



### LOW MAINTENANCE DESIGN

Freddo was designed with low maintenance in mind, with durable materials and easy-to-access components, allowing for less time and effort spent on maintenance tasks. By removing just six screws, it is possible to access all internal parts, ensuring faster and lower-cost repairs, while competitors use modules fixed with rivets, which makes the repair process longer and more expensive.

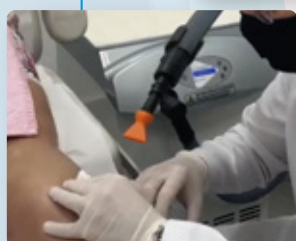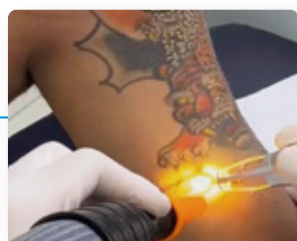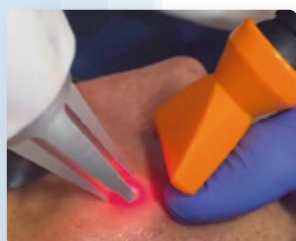

\*Optional articulating arm.





### SUSTAINABLE AND ENVIRONMENTALLY FRIENDLY

All Freddo components are 100% recyclable, contributing to the preservation of the planet and benefiting future generations.



### FREDDO IS ASSEMBLED AND TESTED AT OUR FACTORY IN WINTER PARK, FLORIDA.

You are welcome to contact us and schedule a visit to our factory.



# Freddo

SKIN COOLING DEVICE

TECHNICAL DATA:	
Power supply:	Voltage: 120V-16A ~ 60Hz.
Power consumption:	Max. 11A (1.910 BTU/hour).
Dimensions:	H 28" (71 cm). W 19" (48 cm). L 20" (51 cm).
Weight:	120 lbs (55 kg).
Classification:	Type of protection against electric shock: CLASS I Degree of protection against electric shock: TYPE B
Operating mode:	15:00 – 90:00 min. timer – ADJUSTABLE.



+1 (321) 972-4528

[contact@fabinject.com](mailto:contact@fabinject.com)

7006 Stapoint Court, Suite I, Winter Park, Florida 32792  
Eastpoint Industrial Park.

   /fabinjectllc

[www.fabinject.com](http://www.fabinject.com)